



toto<sup>®</sup>

Single handed  
care, keeping  
your patients  
moving

Automated Patient Turning



Consistent and controlled  
automated lateral patient  
turning system

toto<sup>®</sup>



Time to Care

*for patients, each other  
and yourself*

## Toto Lateral Turning System

### Toto automated lateral turning system

Patients who are immobile are at the highest risk of developing pressure ulcers. Frequent repositioning redistributes pressure and helps prevent pressure damage.

The unique Toto turning system is designed to keep your patients moving by automatically turning your patients at regular intervals day and night.

Keeping patients moving is just one step of the ten-step care plan, aSSKINg<sup>1</sup>.



#### 1 Clinically Effective

Toto has been shown to effectively redistribute pressure away from 'at risk' areas but also increases patient and staff compliance.

#### 2 Promoting Safety

Toto removes the need for additional manual handling aids which can cause damage to other devices (such as mattress covers) and patients skin through friction and shear. Toto can also help reduce the risk of personal injury common to caregivers who regularly handle patients

#### 3 Adaptive Profiling Platform

The new Toto Platform is designed for compatibility with any bed frame in any profiled or 'rise and recline' position, allowing for universal installation of the product.

Toto Platform & Toto Cover



#### 4 Ensuring Compliance

Automated turning ensures that patients are kept moving in busy healthcare environments, removing issues such as; other priorities, missing ineffective reminders, being short of staff and monitoring patient positions.

#### 5 Quick to Install and Easy to Use

Weighing less than 6.3kg, delivered fully assembled in a transport bag and ready for use, the new platform can be transported and installed by only one person with no specialist skills or training required.

Lightweight & Rollable Platform



#### 6 Consistent and Controlled Turning

Toto inflates evenly, smoothly and consistently, fully supporting patients from head to toe avoiding the risk of twisting. Even while patients are sleeping, Toto continues to deliver frequent repositioning.



Adaptive profiling platform



! 250kg User Weight Limit

? Toto lateral turning system has a two year warranty



Toto Platform Transport Bag

## Toto Digital Control Unit


The Toto Touch digital control unit provides customised care, with a user determined side selection feature and preset inflation and rest durations.

Also present, are sophisticated blockage and air leak alarms ensuring patient safety, as well as an interface lock to avoid unwanted setting changes.

A convenient therapy pause mode enables patient interaction without having to stop therapy and a service indicator to maintain performance also comes as standard.



Toto automated turning cycle

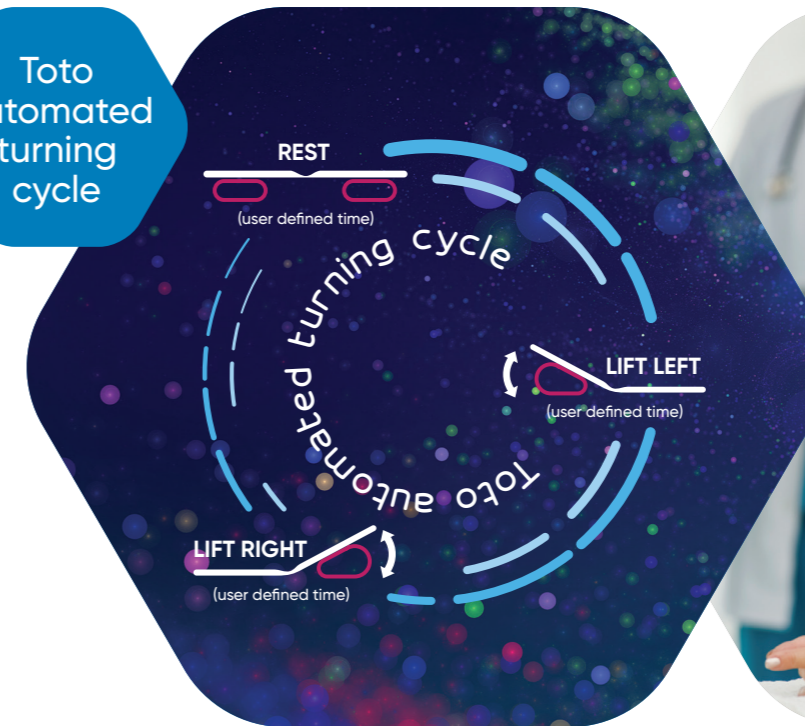

### Automated turning cycle

User determined cycle with independent left and right inflation periods. Five preset time intervals for inflation and rest periods from 30 minutes to four hours, allows for customised care. The system also includes a default mode comprising left and right inflation plus rest period all set at two hours.




### Power on / off



### Therapy pause mode

Therapy pause mode can be engaged by pressing the Power button once. This suspends pump functions enabling patient interaction without having to stop therapy, for example when feeding or undertaking hygiene procedures. The currently selected side selection, inflate and rest time LEDs will flash when activated.

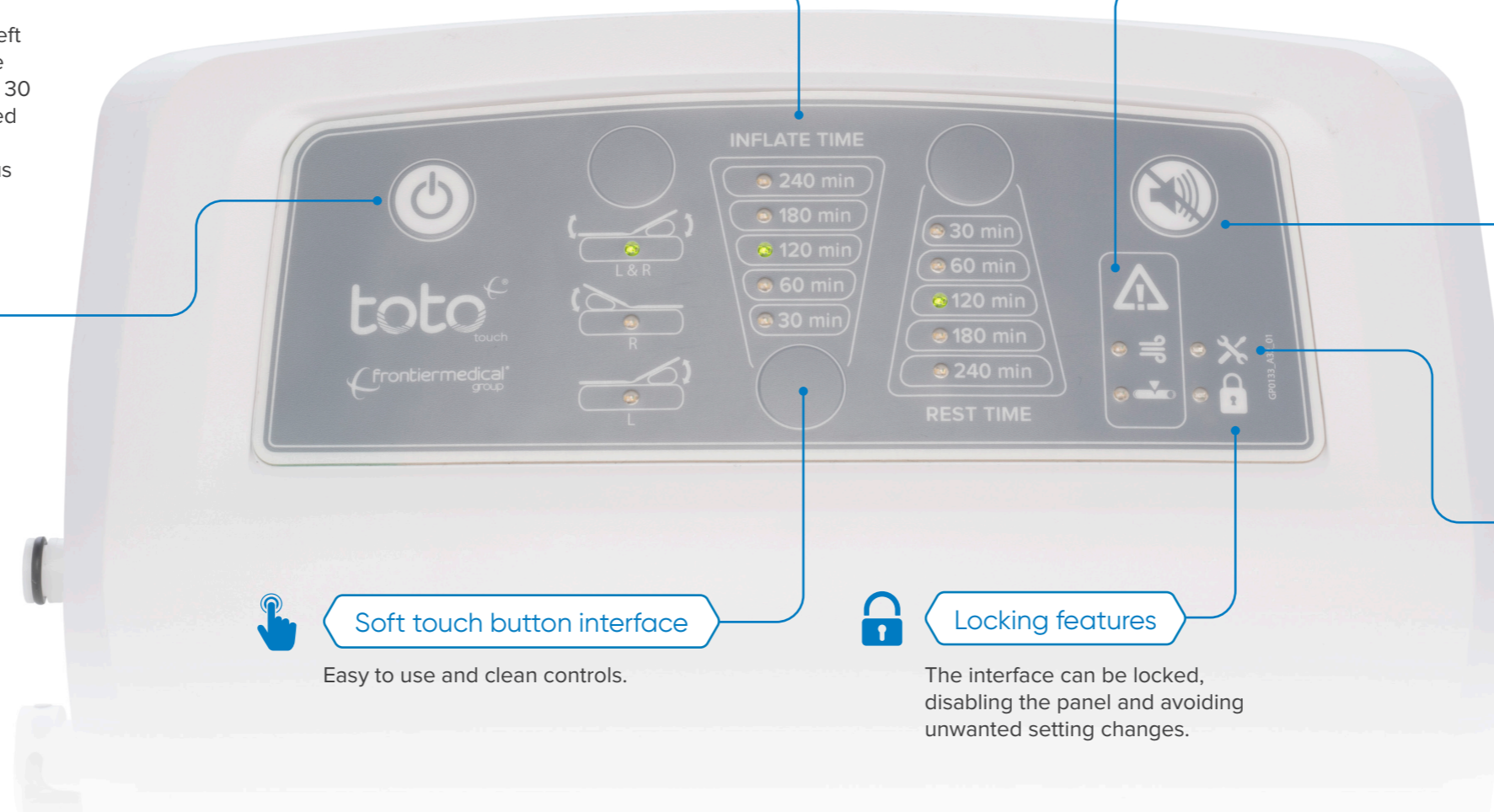

### Soft touch button interface

Easy to use and clean controls.



### Locking features

The interface can be locked, disabling the panel and avoiding unwanted setting changes.

### Alarm features



The control unit incorporates five audible and visual alarms including blockage, low pressure, power loss and service indicators.

### Audio / visual indicators



A series of LEDs and an audible alarm indicate chosen functional options or when trouble shooting needs to be undertaken. A mute function allows for the temporary suspension of an alarm.

### Service indicator



An amber LED indicates when routine servicing is due to ensure optimal performance of the unit.



## Toto Alternative Uses



### Supporting pain management

Toto lateral turning system can be used for those patients who experience pain when manually repositioned.

The gentle tilting action means that the patient can still receive the benefits of repositioning to an angle of approximately 30 degrees without the need for manual intervention which may cause or exacerbate discomfort.



### Repositioning distressed or aggressive patients

For those patients who become distressed or aggressive when repositioned manually, the toto lateral turning system can offer a potential answer to the question of how to position these patients safely.

It will tilt the patient smoothly and gently so that they are less likely to become distressed or aggressive than they would be with manual intervention and thereby ensuring that they still have all the benefits that tilting brings.



### Respiratory management

Moving and tilting patients with respiratory conditions is a common intervention to help manage excess secretions, and assist with postural drainage.

Toto can be used alongside other interventions to aid this process by utilising the tilting action of the system.



### End of life and palliative care

Patients at the end of life are particularly susceptible to pressure ulceration NPIAP/EPUAP/PPPIA 2019, as such, they need to be regularly repositioned to reduce this risk of pressure damage.

This is often upsetting for family and carers and difficult for those looking after them. The Toto lateral patient turning system can be used to ensure that the patient is repositioned and skin management is achieved throughout this difficult time.



## Toto Rental



### Special Products and Service

#### Indications

The Toto lateral turning system is suitable for the prevention of pressure ulcers for the following patients:

- Identified as being at moderate to high risk of developing pressure ulcers
- Weighing up to 250kg (551lbs)
- Who are unable to change their position without assistance
- Identified as requiring regular turning and who are non-compliant with manual turning schedules

#### Contraindications

- Do not use with patients with unstable spinal cord injuries or who are undergoing traction
- Do not use with patients with equilibrium disturbances
- Not to be used where patients must be supported by wedges or pillows etc

#### Considerations

- The use of bedsides and bed rails is recommended
- Assess for the risk of entrapment

#### Storage

- When not in use, remove the cover to prevent damage
- Place a protective covering over the platform and do not store directly on the floor
- Never store other items on top of the Toto platform
- Do not store next to radiators or other heating devices
- Do not allow sharp objects to penetrate the cover

## Order Codes

Description	Dimensions	Order Code
Toto Platform	L 1850mm x W 745mm	4100000
Toto Touch Control Unit	L 240mm x W 100mm x H 200mm	4000000



### Special Products and Service

---

For more information:

Units 5 and 6 Lower Vicarage Road  
Southampton  
Hampshire  
SO19 7RJ  
T: 02380 420678  
E: [info@easycare systems.co.uk](mailto:info@easycare systems.co.uk)  
W: [www.easycare systems.co.uk](http://www.easycare systems.co.uk)

