



MANGAR Camel Lifting Cushion

The Camel Lifting Cushion will both sit up and lift a fallen person.

EASY TO FIT

- Lifting cushion is inflated using an Airflo Compressor with a simple hand control
- Simple to use, requiring very little training
- Suitable for indoor and outdoor use
- Can be used independently or assisted

SIMPLE TO MAINTAIN

- Clean using a proprietary non-abrasive mild liquid cleaner.
- Easy to deflate, fold away and store in the carry bag provided
- Does not need a LOLER certificate
- The material is designed to exceptionally high standards with extraordinary durability, strength & performance in the toughest of environments

SAFE TO USE

- Lifts the frailest person up to 320 kg from the floor
- Stable and supportive
- Ideal for carers of individuals with neurological disorders who fall regularly

IMPROVES WELL-BEING

- Provides a safe, dignified lift for a fallen person
- Reduces the risk of musculoskeletal injury for the carer
- Individuals with a cognitive impairment are less distressed when lifted with a Camel
- The backrest gives additional support and is particularly suitable for bariatric patients

TECHNICAL SPECIFICATIONS

| | |
|-----------------|--------|
| Max user weight | 320 kg |
|-----------------|--------|

| | |
|------------|-------|
| Max height | 56 cm |
|------------|-------|

| | |
|------------|------|
| Min height | Flat |
|------------|------|

| | |
|------------|-------|
| Seat depth | 47 cm |
|------------|-------|

| | |
|--------------|--------|
| Total length | 112 cm |
|--------------|--------|


| | |
|-------|-------|
| Width | 70 cm |
|-------|-------|

| | |
|----------------------|--------|
| Total product weight | 6.5 kg |
|----------------------|--------|

| | |
|----------------------|-------|
| Back/headrest length | 81 cm |
|----------------------|-------|

| | |
|-------------------|----|
| Airflo Compressor | 24 |
|-------------------|----|

| | |
|--------------------------------|---------------|
| Lifting performance per charge | up to 4 lifts |
|--------------------------------|---------------|

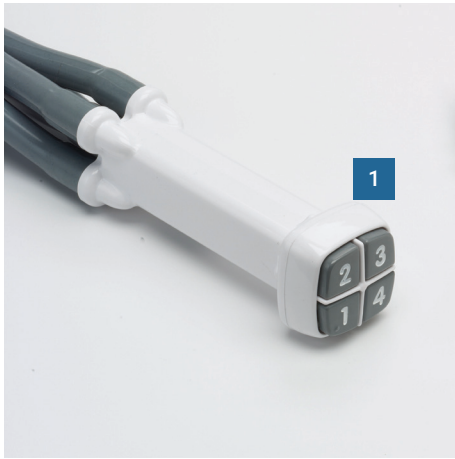

Features

- 1. Easy to setup** - connect the four hoses from the hand control to the corresponding numbers and coloured connectors on the Camel.
- 2. The head should lie on the Camel badge** and the body should lie centrally.
- 3. The backrest** may be partially inflated before the lower sections.
- 4. Made from very strong polyurethane coated nylon material.**
- 5. Each layer inflates individually** ideal for single-handed care.

Accessories

- Slidesheet
- Camel Carry Bag

