



Repose Foot Protector Range



Repose Foot Protectors

Designed specifically to minimise the risk of pressure damage to heels.

repose[®]

Repose® Foot Protector

The Repose Foot Protector and Foot Protector Plus are designed specifically to minimise the risk of pressure damage to heels.

A significant proportion of pressure ulcers occur on the heels due to a combination of pressure and shear and whilst support surfaces assist in reducing pressure,

Smart Valve Technology

The 'Smart valve' ensures Repose products can never be over-inflated following the inflation instructions.


Offloading

Repose Foot Protector and Foot Protector Plus are designed to completely off-load pressure from the heel.

A heel void area incorporated into the Repose Foot Protector provides complete heel offloading.

only a specific device such as Repose Foot Protector and Foot Protector Plus can ensure that pressure on the heel is completely offloaded and so minimising the risk of pressure damage.

Clinical studies have shown the Repose Foot Protectors to be effective in a wide variety of settings.

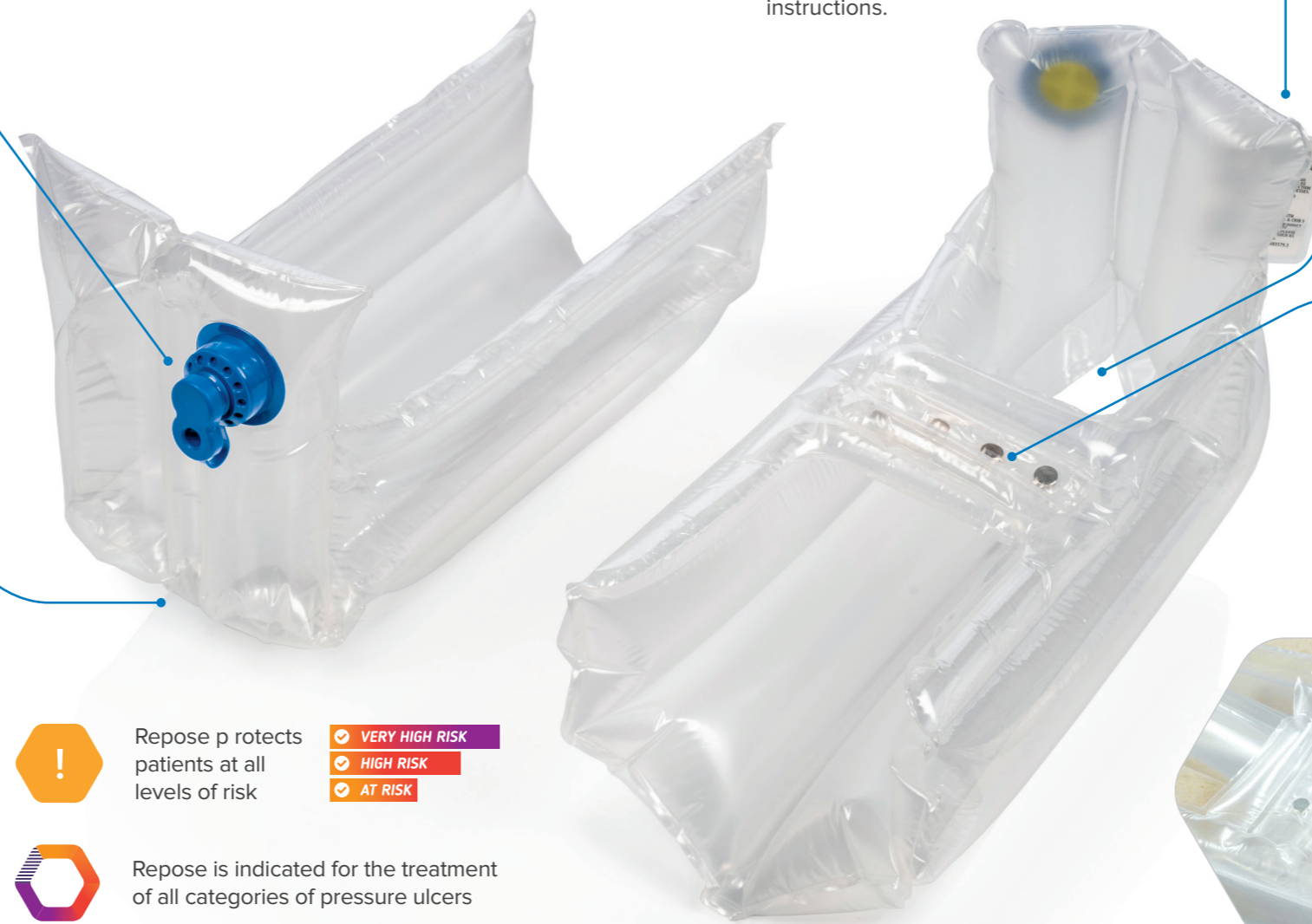

Repose® Foot Protector Plus

The Repose Foot Protector Plus combines the qualities of Repose Foot Protector with the additional benefits of:

ID Label

Provides the option for product traceability plus also includes cleaning and washing instructions.

Ventilation & Visual Inspection Panel

The open heel void of Repose Foot Protector Plus increases ventilation of the foot and so contributes to patient comfort. The open heel void also allows for easier visual inspection of the heel for wound observation and monitoring.


Ventilation & Visual Inspection Panel



Magnaffix Adjustable Strap

Magnaffix is a unique, patented multi-adjustable air filled strap incorporating magnets. This retains the product quickly and easily but is also easy to clean and remove.

The magnets have been sealed within the air cell pocket barrier so that they do not come into contact with the patient's skin to minimise the risk of any pressure or skin tear damage.

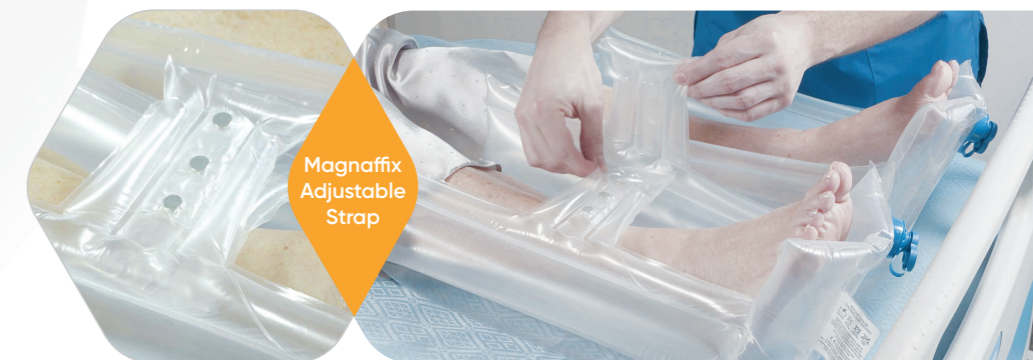

XS Repose Foot Protector Extra Small

Product	1 x Foot Protector Extra Small (Pair) & 1 x Pump
Dimensions	L 290mm x W 200mm x H 180mm

STD Repose Foot Protector Standard

Product	1 x Foot Protector Standard (Pair) & 1 x Pump
Dimensions	L 380mm x W 200mm x H 230mm.

XS Repose Foot Protector Plus Extra Small

Product	1 x Foot Protector Plus Extra Small (Pair) & 1 x Pump
Dimensions	L 290mm x W 200mm x H 180mm

STD Repose Foot Protector Plus Standard

Product	1 x Foot Protector Plus Standard (Pair) & 1 x Pump
Dimensions	L 380mm x W 200mm x H 230mm.

S Repose Foot Protector Small

Product	1 x Foot Protector Small (Pair) and 1 x Pump
Dimensions	L 330mm x W 200mm x H 205mm.

L Repose Foot Protector Large

Product	1 x Foot Protector Large (Pair) & 1 x Pump
Dimensions	L 510mm x W 190mm x H 330mm

S Repose Foot Protector Plus Small

Product	1 x Foot Protector Plus Small (Pair) and 1 x Pump
Dimensions	L 330mm x W 200mm x H 205mm


Sizing Guide

When selecting an appropriate size Repose Foot Protector, you should ensure to measure from the base of the heel to the mid-point of the back of the knee, whilst the patient is sitting.

This is shown in the diagram below. Correspond your measurements to the leg size chart to select an applicable product.

Product Size	Leg Measurement
XS Extra Small	460mm - 410mm
S Small	500mm - 450mm
STD Standard	550mm - 490mm
L Long	690mm - 550mm

Maximum calf circumference 330mm

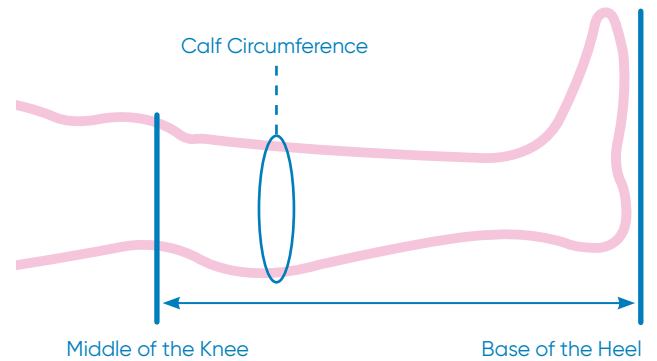

easyCare
systems

Special Products and Service

For more information:
Units 5 and 6 Lower Vicarage Road
Southampton
Hampshire
SO19 7RJ
T: 02380 420678
E: info@easycare systems.co.uk

repose[®]